



omniPAGE Large Format

Ideal for a variety of large format vertical gel applications, these Large Format vertical gel systems offer advanced features for enhancing gel resolution and ease of use, essential when handling gels of this size.



Each unit contains ultra-soft silicone seals for easy plate sealing and trouble-free runs, even over extended run times. Resolution is enhanced by using an aluminium heat sink plate, essential for even sample migration. Added convenience is provided by a removable lower buffer tank and upper buffer drainage tap.

Special buffer chambers allow either low buffer volumes to be used for economy or high buffer volumes to be used for extended runs.

A wide range of interchangeable comb and spacer options allows many techniques to be easily accomplished including; DNA Sequencing, 2-D analysis, Microsatellite analysis, DNA fingerprinting, Gel shift assays, Single-Strand Conformation Polymorphism (SSCP), Heteroduplex and Oligonucleotide analysis.

KEY FEATURES

- Run up to 96 samples
- Enhanced gel heat homogenisation
- Variable low or high buffer volumes
- 20 x 50cm or 33 x 45cm formats

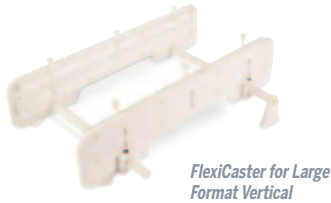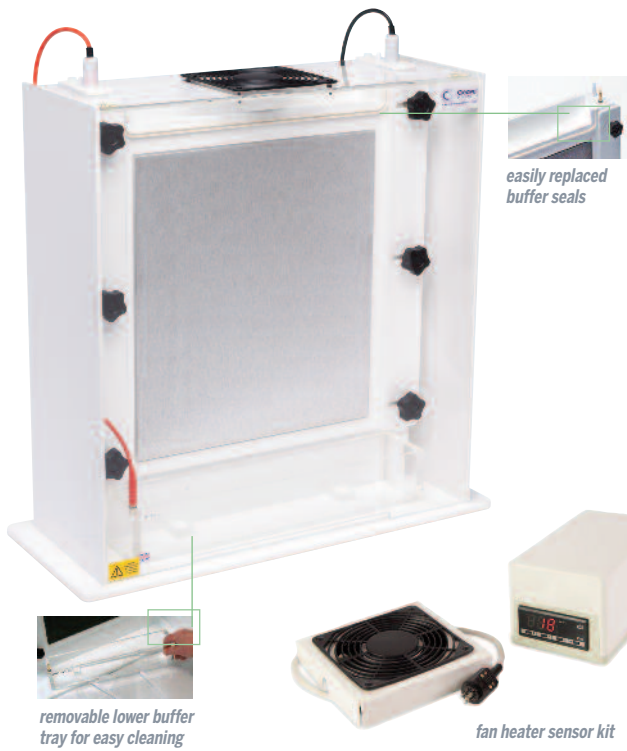

These sturdy racks are designed for safe drying and storage of glass plates. The small rack can hold up to 20x 2mm thick plates while the larger rack can accommodate up to 10x 5mm thick glass plates.

TECHNICAL SPECIFICATIONS	
Plate dimensions (w x l)	CSQ20 20 x 50cm CSQ33 33 x 45cm
Max sample capacity	CSQ20 48 samples CSQ33 96 samples
Buffer volume	CSQ20 Min 500ml, Max 1000ml CSQ33 Min 800ml, Max 2000ml
Combs available No. of teeth Thicknesses	24, 48, 80, 96 0.25, 0.35, 1, 1.5mm

ORDERING INFORMATION

CSQ20	Large Format Vertical , 20cm wide, glass plates, 0.35mm spacers, 48 sample comb		
CSQ20-NG	Glass plates, Notched, pk/2	CSQ20-S0.35	Spacer set 0.35mm
CSQ20-PG	Glass plates, pk/2	CSQ20-S1	Spacer set 1mm
CSQ20-S0.25	Spacer set 0.25mm	CSQ20-S1.5	Spacer set 1.5mm
CSQ33	Large Format Vertical , 33cm wide, glass plates, 0.35mm spacers, 48 sample comb		
CSQ33-NG	Glass plates, Notched, pk/2	CSQ33-S1	Spacer set 1mm
CSQ33-PG	Glass plates, pk/2	CSQ33-S1.5	Spacer set 1.5mm
CSQ33-S0.25	Spacer set 0.25mm	CSL-FHS	Fan heater sensor kit for large format vertical units RRCSQ20 and RRCSQ33
CSQ33-S0.35	Spacer set 0.35mm		
CSL-MGPR	Mini Glass Plate Rack for 20x 2mm Plates	CSL-LGPR	Large Glass Plate Rack for 10x 5mm Plates

CODE	DESCRIPTION	SAMPLE VOLUME PER WELL	CODE	DESCRIPTION	SAMPLE VOLUME PER WELL
CSQ20-0.25-24	Comb 24 sample, 0.25mm thick, Sharks tooth	7µl	CSQ20-1-24SQT	Comb 24 sample, 1mm thick, Square tooth	40µl
CSQ20-0.25-48	Comb 48 sample, 0.25mm thick, Sharks tooth	3µl	CSQ20-1-48SQT	Comb 48 sample, 1mm thick, Square tooth	20µl
CSQ33-0.25-48	Comb 48 sample, 0.25mm thick, Sharks tooth	7µl	CSQ33-1-48SQT	Comb 48 sample, 1mm thick, Square tooth	35µl
CSQ33-0.25-96	Comb 96 sample, 0.25mm thick, Sharks tooth	3µl	CSQ33-1-80SQT	Comb 80 sample, 1mm thick, Square tooth	20µl
CSQ20-0.35-24	Comb 24 sample, 0.35mm thick, Sharks tooth	9µl	CSQ20-1.5-24SQT	Comb 24 sample, 1.5mm thick, Square tooth	60µl
CSQ20-0.35-48	Comb 48 sample, 0.35mm thick, Sharks tooth	5µl	CSQ20-1.5-48SQT	Comb 48 sample, 1.5mm thick, Square tooth	30µl
CSQ33-0.35-48	Comb 48 sample, 0.35mm thick, Sharks tooth	9µl	CSQ33-1.5-48SQT	Comb 48 sample, 1.5mm thick, Square tooth	50µl
CSQ33-0.35-96	Comb 96 sample, 0.35mm thick, Sharks tooth	5µl	CSQ33-1.5-80SQT	Comb 80 sample, 1.5mm thick, Square tooth	30µl